



NTMNEXATM
NPWT System

**INSTRUCTIONS FOR USE
(US VERSION)**

NTMNEXATM
MEDICAL
An AOTI Company

All descriptions and specifications are correct at time of print. Information may be subject to change at any time.

TABLE OF CONTENTS

INTRODUCTION: 4

INTENDED PURPOSE: 4

INDICATIONS AND CONTRAINDICATIONS: 6

PRECAUTIONS: 6

WARNINGS: 8

SETTING UP THE DEVICE: 10

DEVICE AUDIBLE INDICATORS: 11

NEXA™ NPWT SYSTEM AND COMPONENTS: 12

INSERTING THE NEXA™ FLUID CONTAINER PACK: 13

REMOVING THE NEXA™ FLUID CONTAINER PACK: 14

OPERATING THE NEXA™ DEVICE: 15

DAYS REMAINING FEATURE / DEVICE LIFETIME ALERT: 16

USING THE CARRY CASE: 17

APPLYING THE NEXA™ DRESSING 18

REMOVING THE NEXA™ DRESSING 19

TROUBLESHOOTING: 20

CLEANING AND MAINTENANCE / DISPOSAL: 21

TECHNICAL SPECIFICATIONS: 22

EXPLANATION OF SYMBOLS: 23

ELECTROMAGNETIC COMPATIBILITY: 25

ORDERING INFORMATION: 29



WARNING: DO NOT ATTEMPT TO USE THIS DEVICE BEFORE READING THE INSTRUCTIONS FOR USE!

INTRODUCTION:

Description of the NEXA™ NPWT System.

The NEXA NPWT System is an integrated negative pressure wound management system for use in home, acute and extended care settings. The system applies a negative pressure to a sealed wound dressing and promotes wound healing through the removal of exudates into a disposable fluid container.

The system consists of the following key components:

1. NEXA Device: A portable device that contains a pump and a rechargeable Battery and is supplied with a Power Supply, and a Carry Case.
2. NEXA Fluid Container Pack: A single use polymeric flexible fluid container that stores the exudate removed from the wound.
3. NEXA Dressing (-75mmHg or -125mmHg): Sterile single use wound dressing components that are in contact with the wound tissue and a means of sealing to the peri-wound area. Interconnect tubing includes an integrated negative pressure relief valve that limits the maximum pressure able to be applied to the wound and are pre-set and not adjustable by the user.

INTENDED PURPOSE:

The NEXA NPWT system is intended for patients who may benefit from the application of negative pressure to the wound to promote wound healing through the removal of excess exudates, infectious material and tissue debris.

It is intended for use in home care, long-term care and acute settings only when prescribed by a Health Care Professional.

NEXA NPWT may be used in a clinical environment such as a hospital or clinic when treatment is under the supervision of a Health Care Professional.

NEXA NPWT foam dressings should not be placed directly in contact with exposed blood vessels, anastomotic sites, organs or nerves

NEXA NPWT may be used on any part of the body with a wound that is not contraindicated.

INTENDED USER:

NEXA Device: The intended user of the NEXA Device is health care providers/carer/ and/or the patient in a hospital setting and the patients home setting.

NEXA Fluid Container Pack: The intended user of the NEXA Fluid Container Pack is health care providers/carer and/or the patient in a hospital setting and the patients home setting.

NEXA Dressing: The intended user of the NEXA Dressing is a qualified health care professional in a hospital setting and the patients home setting.

INDICATIONS FOR USE:

Appropriate wound types to include:

- Chronic wounds
- Pressure ulcers
- Diabetic foot ulcers
- Venous leg ulcers
- Acute wounds
- Surgical incisional wounds
- Subacute and dehisced wounds
- Partial thickness burns, flaps and grafts

CONTRAINDICATIONS FOR USE:

The device is not recommended for treatment of the following conditions:

- Presence of necrotic tissue
- Malignancy in wound
- Untreated Osteomyelitis
- Untreated malnutrition
- Exposed arteries, veins, nerves, or organs
- Use over anastomotic sites
- Non-enteric and unexplored fistulas
- Necrotic tissue with eschar

PRECAUTIONS:

Precautions should be taken for patients with the following conditions:

- Patients on anticoagulation medicine or who have active bleeding or who have difficult wound haemostasis should be treated with caution. These patients are at an increased risk for bleeding and bleeding complications and should be treated and monitored by properly trained medical caregivers in a controlled setting.
- Wounds that are in close proximity to blood vessels, organs, muscle, and fascia: All blood vessels, organs, muscles, and fascia that are in close proximity to the wound site and/or are exposed and /or are near the skin surface should be properly protected prior to initiating therapy. Patient with infections in the wound and or other parts of the body have to receive proper systemic treatment.
- Irradiated vessels and tissue. These patients are at an increased risk for bleeding and bleeding complications and should be treated and monitored by properly trained medical caregivers in a controlled setting.
- Bony fragments. Sharp edges from bony fragment may puncture blood vessels, organs, muscles, and fascia and may lead to bleeding. Proper care should be taken to cover the bony fragments and protect the wound area and other areas from bleeding.
- Check wound dressing periodically to ensure there is no build up of wound fluid or evidence of any bleeding.
- Dressings should be routinely changed every 48 - 72 hours but no less than 3 times a week or when the wound shows signs of infection.
- Do NOT continue using the system if a Haematoma is observed when changing the dressing.

PRECAUTIONS:

- Spinal Cord Injury: In the event a patient experiences autonomic dysreflexia, discontinue the use of the NEXA™ system and seek advice.
- Bradycardia: The NEXA™ NPWT system must not be placed in proximity to the vagus nerve.
- Enteric Fistulas: The NEXA™ NPWT system is not recommended if enteric fistula effluent management is the sole purpose of management.
- Consider using a skin preparation product to protect the periwound skin. If any signs of irritation appear, cease use of the NEXA™ NPWT System.
- Caregivers should assess patients for adverse reactions on the skin where the device has contact, such as redness (erythema, swelling (edema), irritation, sensitization (delayed Type IV hypersensitivity), allergy, immune response, or other reactions.
- Circumferential Dressing Application: Do not use circumferential dressings unless required to maintain a seal in the presence of anasarca or excessively weeping extremities. It is recommended to use smaller strips to obtain a seal rather than one large drape. Do not stretch the drape upon securement.
- Care must be taken to assess the wound for anastomotic sites. If present, use a non-adherent interface layer between the wound and the foam.
- To prevent tissue ingrowth, consider using an interface layer between the foam and wound. Also consider selecting a lower pressure dressing.
- The NEXA™ NPWT System contains a magnet to facilitate the latch engagement when closing the cover. Caution should be taken when using this device within 6 inches of other magnetically susceptible devices.

GENERAL PRECAUTIONS FOR ALL INDICATIONS FOR USE:

- It is important that a doctor, nurse or other qualified healthcare provider evaluates the patient to ensure that the use of this device is an appropriate therapy.
- Prior to issuing the system, the patient should be assessed for their knowledge, skill level and amount of training required enabling the patient to use the system themselves.
- To reduce the risk of transmission of blood-borne pathogens, regardless of the diagnosis or presumed infection status, all users should take suitable precautions as defined in locally applicable standard operating procedures for infection control.

WARNINGS:

The following Warning statements describe the potential for serious consequences to the patient, such as death, injury, or adverse reactions. Failure to read and follow all instructions in this manual prior to use may result in death or injury of the patient.

- Do not use any Power Supply to re-charge the device other than that supplied.
- The socket in which the Power Supply is connected must be accessible at all times.
- Only use dressings and accessories that are approved for use with the NEXA™ Device.
- Do not use dressings or accessories that are past expiration date.
- Physician should consider the patient's size and weight when prescribing this device.
- The device is not safe for use with an MRI and must be disconnected from the patient prior to MRI.
- Do not use the device in a Hyperbaric Chamber or in the presence of flammable gases. The patient's dressing may remain in place when disconnected from the unit.
- To prevent unintentional gauze/foam retention, all dressings should be carefully removed from the wound and the entire wound bed. Upon removal of the dressings, the wound bed should be cleaned in accordance with standard wound care practices (or facility guidelines), prior to the application of new sterile dressing.
- Ensure that there are no pockets left in the wound or wounds after application of the dressings.
- Infected wounds must be inspected more frequently for signs of increased infection or sepsis.
- Patients who do not have adequate haemostasis, and who are currently taking anticoagulation or platelet aggregation inhibitor therapy, have an increased risk of bleeding with or without the NPWT device.
- Defibrillation: Remove the NEXA™ Dressing if defibrillation is required in the area of dressing placement. The dressing may inhibit transmission of electrical energy and/or patient resuscitation.
- All arteries, veins, tendons, ligaments, nerves, and organs must be covered completely prior to application of the device.
- Use with caution on any patient with increased risk of bleeding due to the presence of weakened or friable blood vessels or organs, conditions when suture of the blood vessels has taken place or when there is (localised to the wound) any infection, trauma, or radiation, as if not controlled well, could potentially be fatal.
- Infected tissue such as blood vessels may have a weakened structure and have to be treated with care. Infected blood vessels may bleed more readily than normal blood vessels.
- There is a risk of strangulation or asphyxiation from tubing, cables and the Carry Case strap. Retain excess tubing to patients body using supplied drape.

WARNINGS:

- In the event of a device malfunction return the equipment to the authorised service centre.
- There are no serviceable parts inside the therapy unit therefore do not open the device.
- No modification of the device is allowed!
- The device must never be used to remove explosive gases and flammable or corrosive fluids.
- The device must not be operated in damp rooms or when taking a bath or shower. Avoid moisture on plug and switches.
- Never plunge the device into water or liquids, not even when it is switched off.
- The unit must not be operated in splash water range, near sources of steam or in locations where there is a danger of explosion.
- EMC Statement: Although the device is compliant with the current EMC regulations applicable, the device may be susceptible to EMC radiation from medical emitters such as electrocautery, MRI, electrosurgery and diathermy devices and action should be taken to avoid exposure.
- The NEXA™ NPWT System is a medical device, not a toy. Keep away from children, pets and pests as they can damage the NEXA™ Device, NEXA™ Dressing and Fluid Container Pack and potentially affect their performance. Keep the system free from lint and dust as they may cause damage and affect performance.
- Keep the system away from sources of heat such as heaters, fireplaces or direct sunlight.

- To ensure the system is working correctly, the dressing and device should be checked regularly throughout the day to check for fluid mobility in the tubing and that the green led is lit.

- The NEXA™ device does not have audible or visual alarms to alert users to a leak or blockage in the system. Please refer to the troubleshooting section to determine how a leak or blockage may be detected and resolved.

- The foam and drape dressing components, that come into patient contact, may cause an allergic reaction. Ensure the foam does not overlap on to healthy skin and always monitor the skin for any signs of reaction.

- Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the member state in which the user and/or patient is established.

SETTING UP THE DEVICE:

- Do not attempt to use the device until reading the instructions.
- The device should only be administered by persons who have been trained in its operation according to the instruction guidelines issued by the supplier or qualified medical staff.
- Before using the device, please read the indications, warnings, precautions and contraindications.
- Check function of the unit prior to use.
- Fully charge the NEXA™ device prior to first use. charging takes approximately 4 hours.
- The device has a sleep mode. If the charging light is not flashing, turn the device on and off again to wake it.
- For optimum performance and safety, charge the device in an ambient temperature of 20°C +/-5°C (68°F +/-9°F)
- Prior to connecting to the power supply adapter, check whether the voltage corresponds with the in building voltage.
- Never connect the Power Supply to defective power sockets. Keep power supply adapter and cable away from external heat sources.
- The device should not be charged or started up if:
 - The power cord or plug are defective
 - The device is not functioning properly
 - The device has been damaged / dropped
 - The device has been dropped into water
 - Obvious defects might restrict safe operation

In any case, remove the Power Supply from the electrical socket and have the unit checked by qualified personnel authorised by AOTI Ltd.

- The device must be used within the Carry Case at all times except to replace the Fluid Container Pack.
- It is the responsibility of the clinician or trained caregiver to determine if the patient's condition allows for mobile use.
- The device can be charged whilst in its Carry Case by pulling back the tab near the socket and inserting the plug.
- Operation of the device is possible while the battery is charging.
- Pay attention to the ambient conditions described in the technical data.
- If the device is operated at ambient temperatures outside the stated temperature range (see "Technical Data"), the performance may be reduced and the unit or the electronics and battery may get damaged.
- If the device has been stored for a pro-longed period without use then the battery should be checked by re-charging fully in accordance with the Instructions, and if necessary returned to a qualified organisation to replace the battery.


DEVICE AUDIBLE INDICATORS:

The following statements describe the audible indicators that may be observed in the following conditions during device operation.


- **Device Power On**
When the Power button is pressed to turn the device ON the audible sounder will beep once.
- **Fluid Container Full**
When the fluid container full pressure switch is activated, the audible sounder will beep for 10 seconds duration (pulsed beep), repeating every 15 minutes.
- **Fluid Container Cover Open**
When the fluid container cover is opened during operation, the audible sounder will beep for 10 seconds duration (pulsed beep), repeating every 15 minutes.
- **Battery Low Indicator**
When the battery low indicator is activated, the audible sounder will beep for 5 seconds duration (pulsed beep), repeating every 45 minutes.
- **Device Lifetime Alert**
When the device end of life indicator is activated, the audible sounder will beep for 5 seconds duration (continuous beep), repeating every 24 hours.

NEXA™ NPWT SYSTEM AND COMPONENTS:

NEXA™ Device

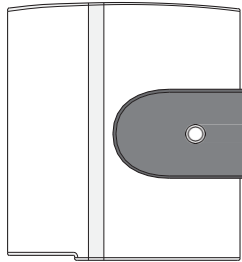

NEXA™ Power Supply

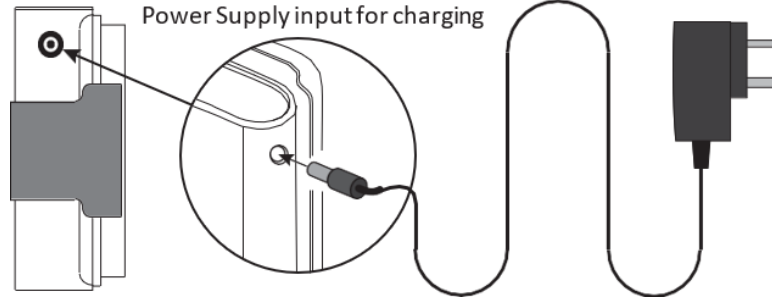


NEXA™ Carry Case




NOTE: THE DEVICE HAS A SLEEP MODE. IF THE CHARGING LIGHT DOES NOT FLASH WHEN CONNECTED TO THE MAINS POWER, TURN THE DEVICE ON AND OFF AGAIN TO WAKE THE DEVICE FROM ITS SLEEP MODE.

Consumables (sold separately)

NEXA™ Fluid Container/tubing/Connector

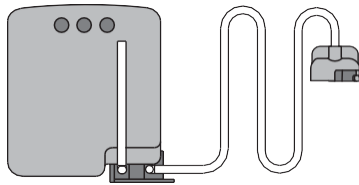



NEXA™ Connector/tubing/pad, Foam, Drape (-75mmHg or -125mmHg)

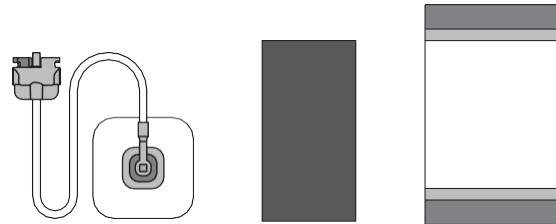



INSERTING THE NEXA™ FLUID CONTAINER PACK:

- NOTE:**
- IT IS RECOMMENDED TO REPLACE THE FLUID CONTAINER TWICE A WEEK TO ENSURE OPTIMUM THERAPY LEVELS AND REDUCE ODOUR.
 - IN ANY CASE IT MUST BE CHANGED ONCE A WEEK, IRRESPECTIVE OF FLUID LEVELS.



OPEN

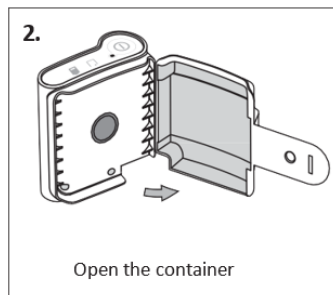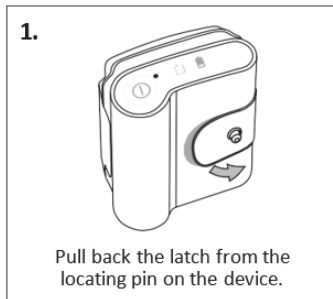

INSERT

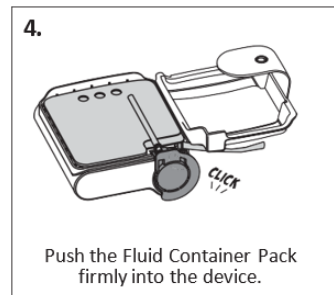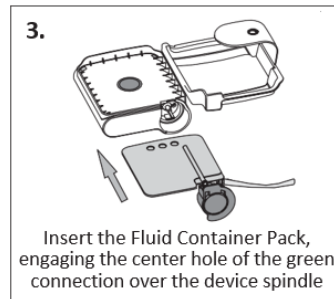



CLOSE

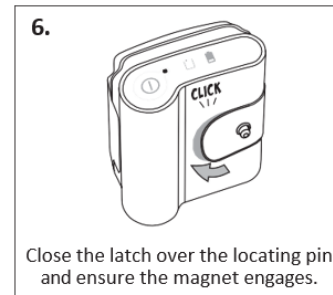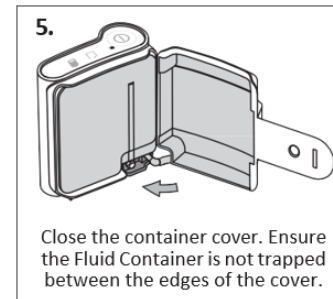



REMOVING THE NEXA™ FLUID CONTAINER PACK:

NOTE: DISPOSE OF THE USED FLUID CONTAINER PACK ACCORDING TO INSTITUTION AND LOCAL ENVIRONMENT REGULATIONS.



TURN OFF/DISCONNECT

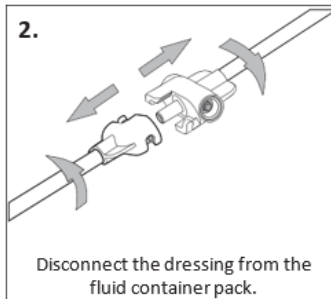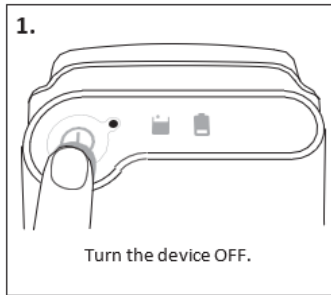

OPEN

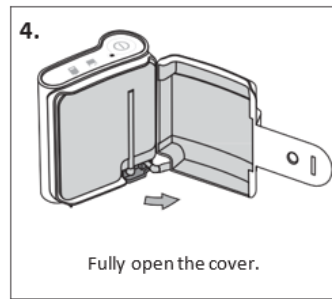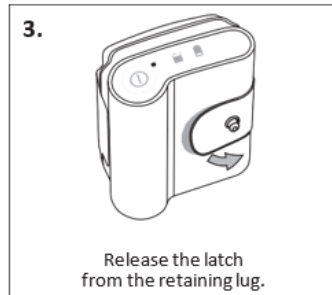


RELEASE/REMOVE

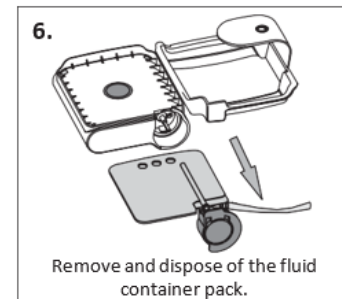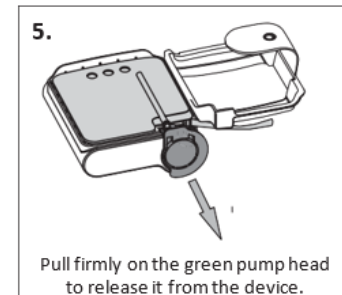

OPERATING THE NEXA™ DEVICE:

NOTE: DO NOT TURN THE DEVICE ON WITHOUT A FLUID CONTAINER PACK INSERTED AND A NEXA™ DRESSING CONNECTED

TURNING ON THE DEVICE

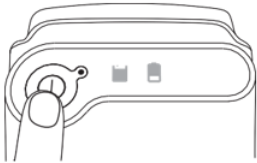

Press the Power button to turn the device ON. The green LED will illuminate, and the audible sounder will beep once.

If the green LED does not illuminate, charge the device immediately using the Power Supply provided and turn the device ON to wake it from sleep mode. The orange light will flash to indicate charging. Charging takes approximately 4 hours. The device will run for approximately 10 hours with a fully charged battery.

The device will start at maximum speed to help apply a negative pressure to the dressing quickly. Observe the dressing and foam to ensure there is a good seal and the foam compresses.

After 3 minutes, the device AUTOMATICALLY reduces its speed. This will be the normal operating speed.

To turn the device OFF, press the power button once and the green LED will turn off.

LED LIGHT ALERTS AND THE USER INTERFACE



Green LED is lit = **Power is ON**

Therapy is ACTIVE



Green LED is lit = **Power is ON**

Orange LED is lit = **Battery LOW**

Therapy is ACTIVE



Green LED is lit = **Power is ON**

Orange LED flashes = **Battery CHARGING**

Therapy is ACTIVE



Green LED is lit = **Power is ON**

Red LED is lit = **Container is FULL or Latch is not engaged**

Therapy is INACTIVE


DAYS REMAINING FEATURE:

The device has a feature to see how many days remain on the device. Open the grey cover. Press and hold the grey button and at the same time press and hold the power button for 8 seconds. Do not release the buttons until the LEDs start flashing. The RED light will flash to indicate the 10-day increments and the ORANGE light will flash to indicate the 1-d day increments.

For example:


- 4 Red flashes and 2 Orange flashes will indicate 42 days remaining.
- 2 Red flashes and 8 Orange flashes will indicate 28 days remaining.

To resume therapy, press the power button.

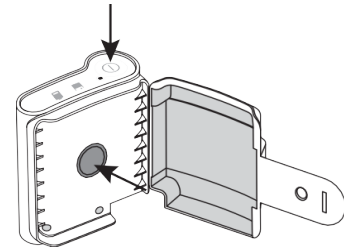


DEVICE LIFETIME ALERT:

The device has a lifetime of 1344 hours (8 weeks) continuous use. The device will alert the user when days are remaining.

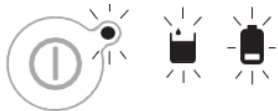

At the beginning of the 3-day warning period, all LED's will flash for a short period of time. The audible sounder will beep for 5 second duration and repeat every 24 hours.

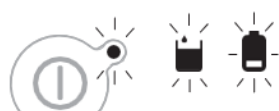

Every hour all LED's will flash for a short period of time.



At the end of device lifetime (1344 hours), all LED's will flash for a short period of time and the device will shut down.




If the device is switched on following shut down, all LED's will stay illuminated but the device will not function.

If the device is switched off and on again during the 3 day warning period, all LED's will flash for a short period of time.


USING THE CARRY CASE:

The carry case can be attached to the patient using the shoulder strap or belt strap.
To charge the device, push the Power Supply connector through the opening in the carry case into the socket of the device.
Place on a table to remove the device from the device from the carry case.

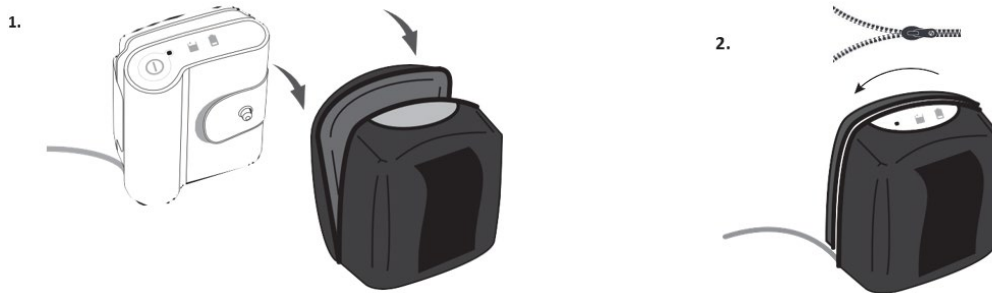









POSITIONING THE NEXA™ NPWT SYSTEM:

Keep the device in the upright position when it is operating.
This will enable optimum performance and allow for the maximum 120ml of fluid to be contained in the Fluid Container.

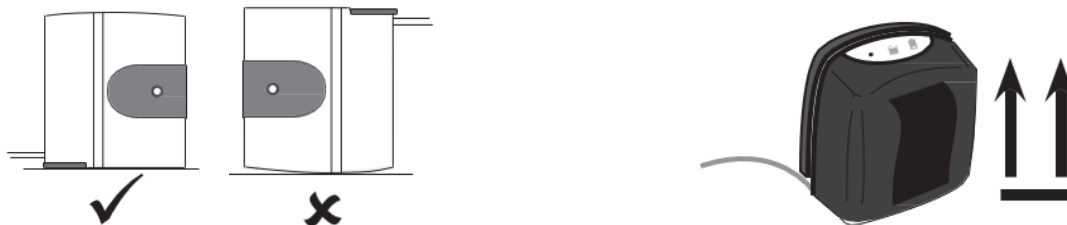

APPLYING THE NEXA™ DRESSING (TO BE PERFORMED BY A CLINICIAN OR HEALTHCARE PROFESSIONAL ONLY):

- NOTE:**
- ALWAYS ENSURE THE DEVICE IS SWITCHED OFF WHEN CONNECTING THE DRESSING.
 - SELECT THE APPROPRIATE DRESSING PRESSURE VALUE, -75mmHg OR -125mmHg, AS CLEARLY LABELLED ON THE PACKAGING.
 - DO NOT USE THE STERILE DRESSING IF THE PACKAGING HAS BEEN DAMAGED OR IS ALREADY OPEN.
 - DRESSINGS SHOULD BE CHANGED EVERY 48 - 72 HOURS BUT NO LESS THAN 3 TIMES A WEEK. NOT TO BE RE-STERILISED OR RE-USED.

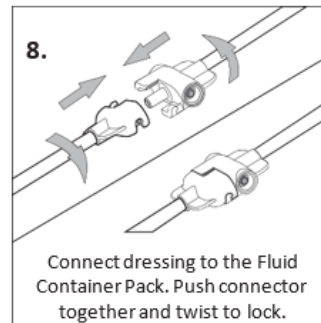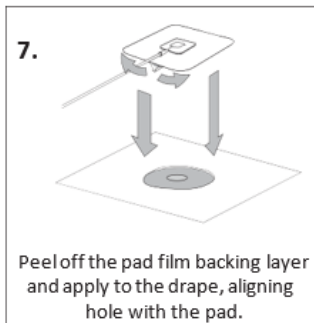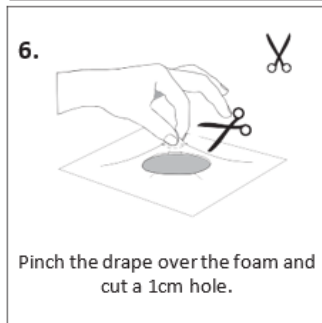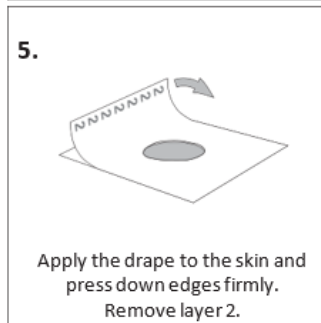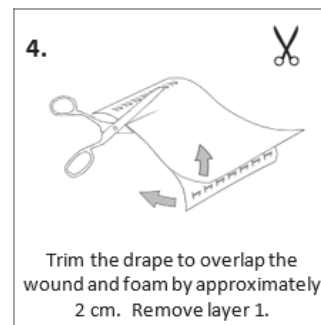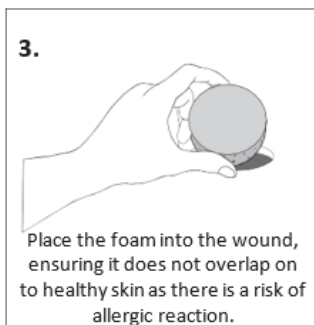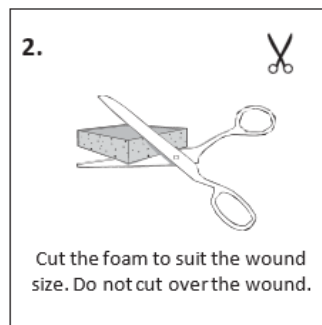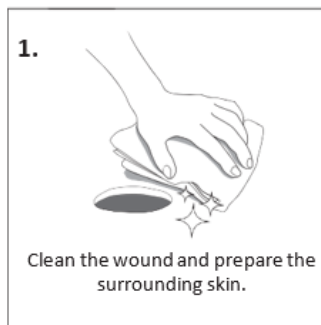

REMOVING THE NEXA™ DRESSING (TO BE PERFORMED BY A CLINICIAN OR HEALTHCARE PROFESSIONAL ONLY):

NOTE: SOME PATIENTS MAY EXPERIENCE PAIN UPON DRESSING REMOVAL. REMOVE DRAPES GENTLY.

1. TURN THE NEXA™ DEVICE OFF.
2. DISCONNECT THE DRESSING FROM THE FLUID CONTAINER PACK.
3. GENTLY STRETCH THE DRAPE HORIZONTALLY TO RELEASE THE ADHESIVE DRAPE FROM THE SKIN.
4. DO NOT PULL THE DRAPE VERTICALLY AS THIS WILL CAUSE MORE PAIN.
5. GENTLY REMOVE THE ENTIRE DRAPE FROM THE PATIENT.
6. GENTLY REMOVE THE FOAM FROM THE WOUND ENSURING THAT NO SMALL FOAM PIECES ARE LEFT IN THE WOUND.
7. CLEAN RESIDUE REMAINING FROM THE DRESSING COMPONENTS IN THE PERIWOUND AREA.
8. DISPOSE OF THE USED DRESSING ACCORDING TO INSTITUTION AND LOCAL ENVIRONMENTAL REGULATIONS.

BEFORE APPLYING A NEW DRESSING ENSURE THE SKIN IS DRY AND CLEAN AND NO FOAM REMAINS IN THE WOUND. FOLLOW THE 'APPLYING THE NEXA™ DRESSING' INSTRUCTIONS.

TROUBLESHOOTING:

LEAKS: There is no leak alarm as the device will continue to remove fluid, even with reduced vacuum, thereby reducing the potential for skin maceration. If the device is unable to pull down the dressing or the dressing is not compressed, check the following:

- Dressing is firmly sealed around the wound. Use additional film drape if necessary around the drape edges.
- Check dressing connector is properly locked.
- Fluid Container Pack is fully engaged into the device.

BLOCKAGES:

If there is no fluid mobility in the tubing and the dressing is not compressed, there may be a blockage. A blockage will also cause the device to make an uncharacteristic slapping noise.

- To resolve a blockage, check that there is no kinking in the tubing or the tube is occluded. If there is, either straighten the kink or replace the dressing and/or Fluid Container Pack.

BATTERY LOW INDICATOR - When the battery reaches low charge, the orange Battery Low light will be lit and the sounder will beep to indicate approximately 2 hours of battery life remaining. Connect the power supply to the device and plug into a power outlet. Once charging, the light will flash intermittently to indicate charging is occurring. Charge time is approximately 4 hours. The device can be used whilst charging. When the device is fully charged, the orange light will go out. The device will operate on battery only for approximately 10 hours.

FLUID CONTAINER FULL/LATCH ENGAGED INDICATOR

When the red light is lit, the device will stop running.

If the red light is lit and the audible sounder is activated there are two possible reasons:

- The Fluid Container is full and requires changing.
- The latch has not been correctly attached to the device. Ensure the latch is engaged over the fastener on the front of the device.
- If the Fluid Container is full of air, ensure the device is operated in its upright position.
- When a new Fluid Container Pack has been installed and the latch engaged, the red light will go out and the device will restart.
- The device will restart at maximum speed for 3 minutes. If required, open and close latch to give an additional 3 minutes of high speed.

LIFETIME ALERT - If the device shows all the LED's flashing or all LED's permanently on, refer to page 16 for the Device Lifetime Alert.

CLEANING AND MAINTENANCE:

Always unplug the power unit from electrical outlet before cleaning. Ensure any local or institutional regulations on hygiene are complied with. The NEXA™ Device is non-serviceable.

The NEXA™ Device is suitable for use by multiple patients. It is necessary to clean and disinfect the device between patients.

Wear suitable gloves for cleaning / disinfection. Routine cleaning of device can be done by wiping down with damp cloth using disinfectant and water or non-aggressive cleaning material.

WARNING! TAKE SPECIAL PRECAUTION TO ENSURE THAT NO CLEANSING SOLUTION IS ABLE TO PENETRATE INTO THE EQUIPMENT. DO NOT RE-USE THE FLUID CONTAINER PACK. DISPOSE OF PROPERLY ACCORDING TO LOCAL AND INSTITUTIONAL GUIDELINES.

DISPOSAL:

The device is made from various electronics, plastics and a Lithium-Ion battery. When the pump is ready for disposal, facilities should follow the local governing guidelines regarding sanitation of disposed device components.

The used Fluid Container Packs, tubes and dressings should be disposed according to the local or facility guidelines for handling infected or bio-hazardous materials.

None of the items should be disposed together with household or facility refuse. Incorrect disposal can have harmful effects on the environment and public health.

CAUTION! PAY ATTENTION TO COUNTRY-SPECIFIC, LOCAL, AND FACILITY REGULATIONS WITH RESPECT TO DISPOSAL, ESPECIALLY WITH REGARD TO DISPOSAL OF USED BATTERIES.

EXPLANATION OF SYMBOLS:



Refer to Instructions for Use



Manufacturer



Date of manufacture



Expiry Date



MR Unsafe



Do not re-use



Do Not Re-sterilize



Caution



Keep Dry



Fragile



Do not use if package is open or damaged



Keep out of direct sunlight



Humidity Limitations



Temperature Limitations

EXPLANATION OF SYMBOLS:



Single Sterile Barrier System

Rx Only

Prescription only

IP02

Ingress Protection Rating:

No protection ingress of solid objects; Protected against falling drops of water, if the case is disposed up to 15° from vertical.



Reference Number



LOT Number



Medical Device



Method of Sterilisation - Ethylene Oxide



WEEE Directive 2012/19/EU

IP20

Ingress Protection Rating:

Protection against solid objects over 12.5mm; no protection against ingress of liquids.



Unique Device Identifier



Serial Number


ELECTROMAGNETIC COMPATIBILITY:

Guidance and manufacturer's declaration - electromagnetic emissions		
The NEXA™ Device is intended for use in the electromagnetic environment specified below. The customer or the end user of the NEXA™ Device should assure that it is used in such environment.		
Emissions Test	Compliance	Electromagnetic environment - guidelines
RF emissions - CISPR 11 (Radiated and Conducted)	Group 1	The NEXA™ Device is suitable in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
RF emissions - CISPR 11 (Radiated and Conducted)	Class B	
Harmonic emissions EN 61000-3-2	Class A	
Voltage Fluctuations / Flicker emissions EN 61000-3-3	Complies	

ELECTROMAGNETIC COMPATIBILITY:

Guidance and manufacturer's declaration - electromagnetic immunity			
The NEXA™ Device is intended for use in the electromagnetic environment specified below. The customer or the end user of the NEXA™ Device should assure that it is used in such environment.			
Immunity Test	EN/IEC Test Level	Compliance Level	Electromagnetic environment - guidelines
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±15kV air	±8 kV contact ±15 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient / burst IEC 61000-4-4	±2 kV for power supply lines ±1 kV for input/output lines	±2 kV for power supply lines ±1 kV for input/output lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	±1 kV differential mode (line to line) ±2 kV common mode (line to earth)	±1 kV differential mode (line to line) ±2 kV common mode (line to earth)	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines EN 61000-4-11	<5% Ut (>95% dip in Ut) for 0.5 cycle 40% Ut (60% dip in Ut) for 5 cycles 70% Ut (30% dip in Ut) for 25 cycles <5% Ut (>95% dip in Ut) for 5 secs	<5% Ut (>95% dip in Ut) for 0.5 cycle 60% dip in Ut for 5 cycles 30% dip in Ut for 25 cycles >95% dip in Ut for 5 secs	Product has internal battery backup. Mains power quality should be that of a typical commercial or hospital environment. If the user of the NEXA™ Device requires continued operation during power mains interruptions, it is recommended the device be powered from an uninterruptible power supply or a battery.
Power Frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
Note: Ut is the mains voltage prior to application of the test level			

ELECTROMAGNETIC COMPATIBILITY:

Guidance and manufacturer's declaration - electromagnetic immunity			
The NEXA™ Device is intended for use in the electromagnetic environment specified below. The customer or the end user of the NEXA™ Device should assure that it is used in such environment.			
Immunity Test	EN/IEC Test Level	Compliance Level	Electromagnetic environment - guidelines
<p>Conducted RF IEC 61000-4-6</p> <p>Radiated RF IEC 61000-4-3</p>	<p>3Vrms 150kHz to 80 MHz</p> <p>3V/m 80 MHz to 2.5 GHz</p>	<p>3Vrms 150kHz to 80 MHz</p> <p>3V/m 80 MHz to 2.5 GHz</p>	<p>Portable and mobile RF communications equipment should be used no closer to any part of the NEXA™ Device, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.</p> <p>Recommended separation distance $d=1.2VP$ $d=1.2VP$ (80 MHz to 800 MHz) $d=2.3VP$ (800 MHz to 2.5 GHz) where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m).</p> <p>Fields strengths from fixed RF transmitters, as determined by an electromagnetic site survey^a, should be less than the compliance level in each frequency range.^b</p> <p>Interference may occur in the vicinity of equipment marked with the following symbol: </p>
NOTE 1: At 80MHz and 800MHz, the higher frequency range applies.			
NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from objects, structures and people.			
<p>a. Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic survey should be considered. If the measured field strength in the location in which NEXA™ is used exceeds 3V/m, NEXA™ should be observed to verify normal operation. If abnormal operation is observed, additional measures may be necessary, such as reorienting or relocating NEXA™.</p> <p>b. Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m</p>			

ELECTROMAGNETIC COMPATIBILITY:

Recommended separation distanced between portable and mobile RF communications equipment and the NEXA™ Device			
The NEXA™ Device is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the NEXA™ Device can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the NEXA™ Device as recommended below, according to maximum output power of the communications equipment.			
Rated maximum output power of transmitter (W)	Separation distance according to frequency of transmitter in meters (m)		
	150 kHz to 80 MHz d=1.2VP	80 MHz to 800 MHz d=1.2VP	800 MHz to 2.5 GHz d=2.3VP
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum power rating of the transmitter in watts (W) according to the transmitter manufacturer.

Note 1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

ORDERING INFORMATION:

NEXA™ Device Kit 414-00-300-XXX
(XXX indicates the country order code, contact your distributor for your order code)

CONSUMABLES:

NEXA™ Dressing -125mmHg 414-00-501 (Box of 5)
NEXA™ Dressing -75mmHg 414-00-511 (Box of 5)
NEXA™ Fluid Container 414-00-502 (Box of 5)
NEXA™ Drape 414-00-600 (Box of 10)
NEXA™ Foam 414-00-700 (Box of 5)

SPARES:

NEXA™ Carry Case 414-00-400
NEXA™ Power Supply 414-00-450 (US)
with country specific adaptor 414-00-451 (UK)
414-00-452 (EU)
414-00-453 (AU)
414-00-454 (CH)

CONTACT INFORMATION:

Notify your nursing contact if any problems or issues are encountered with the system.
For questions regarding this product, supplies, maintenance and information, contact:
Website: www.aotinc.net Email: customerservice@aotinc.net



AOTI Ltd.
Unit 20 Glenrock Business Park,
Ballybane, Galway, Ireland.

All trademarks designated herein are property of AOTI Ltd and its licensors.
The NEXA™ products are protected by global patents, granted and pending. For patents see www.aotinc.net
All rights reserved.

414-06-002 Rev AC

This page is left intentionally blank

